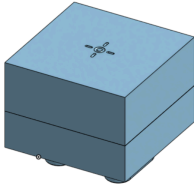


CubXL+ - Calibration Block, *Mounting Instructions*

Revision History.


Rev. No	Title	Image	BOM
0	<i>CubXL+ - Calibration Block</i>		See Appendix 1

3D Printing CubXL+ - Calibration Block

NOTE: The CubXL+ - Calibration Block was originally printed using a *Bambu Labs X1-Carbon* 3D Printer, with *Bambu Lab PLA Basic* filament.

NOTE: The default print settings were used with the exception of support generation being turned ON in order to support overhangs. See Figure 1 for a screenshot of support settings.

Figure 1.

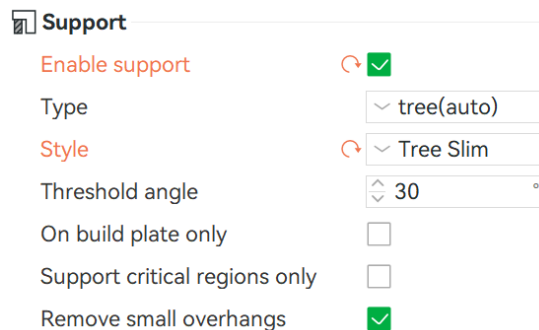

NOTE: The values in orange were changed from the default settings

Installation of CubXL+ - Calibration Block

Required Components.

Item No.	Name	Rev.	Picture	Quantity
1	<i>Calibration Block</i>	0	See Appendix 1	1
2	<i>Calibration Base (CubXL+)</i>	0		1



STEP 1. Align the Calibration Base (CubXL+) in the desired location on the gantry baseplate, ensuring that the two (2) mounting pegs are aligned with two (2) holes on the baseplate of the gantry.

NOTE: Placement of the Calibration Base (CubXL+) is suggested to be in any of the four (4) corners of the gantry baseplate.

STEP 2. Align the four (4) diamond shaped extrusions on the bottom of the Calibration Block with the corresponding four (4) diamond shaped holes on the top of the Calibration Base (CubXL+).

Appendix 1

CubXL+ - Calibration Block_REV. 0 BOM

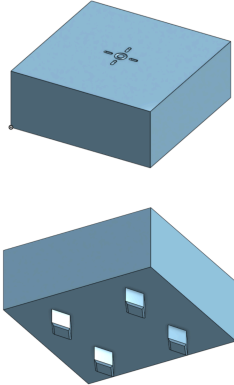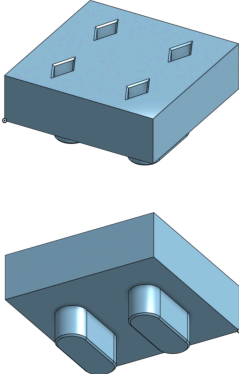Item No.	Title	Image	Revision
1	<i>Calibration Block</i>		0
2	<i>Calibration Base (CubXL+)</i>		0