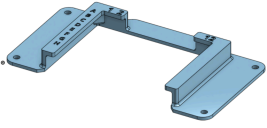
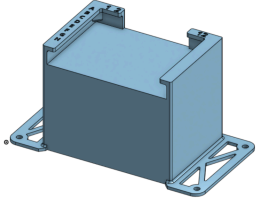
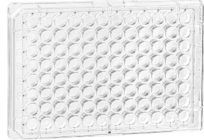


CubXL - ASMI Wellplate Holder, *Mounting Instructions*

Revision History.

Rev. No	Title	Image	BOM
0	<i>CubXL - ASMI Wellplate Holder</i>		N/A
1	<i>CubXL - ASMI Wellplate Holder</i>		N/A

NOTE: The CubXL - ASMI Wellplate Holder was designed to be used in conjunction with the following component(s):


Item No.	Title	Qty.	Image	Link
1	<i>Cell Culture Microplate, 96 Well, PS, F-Bottom</i>	1		https://shop.gbo.com/en/usa/products/bioscience/cell-culture-products/cellstar-cell-culture-microplates/96-well-cell-culture-microplates-clear/655160.html

3D Printing CubXL - ASMI Wellplate Holder

NOTE: The CubXL - ASMI Wellplate Holder was originally printed using a *Bambu Labs X1-Carbon* 3D Printer, with *Bambu Lab PLA Basic* filament.

NOTE: The default print settings were used with the exception of support generation being turned ON in order to support overhangs. See Figure 1 for a screenshot of support settings.

Figure 1.

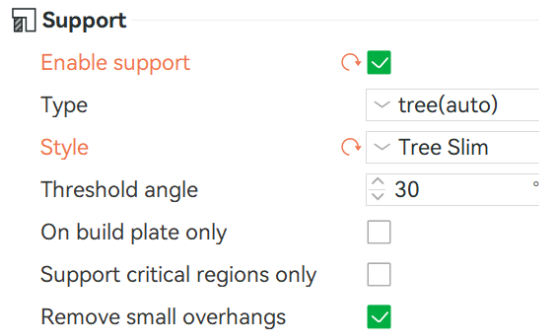

NOTE: The values in orange were changed from the default settings

Installation of CubXL - ASMI Wellplate Holder


Required Components.

Item No.	Name	Rev.	Picture	Quantity
1	<i>CubXL - ASMI Wellplate Holder</i>	0-1	See Revision History	1
2	<i>18-8 SS Button Head Screw, M5, 15mm</i>	N/A		4
3	<i>Zinc Plated Steel Hex Nut, M5</i>	N/A		4

STEP 1. Align the CubXL - ASMI Wellplate Mount in the desired location on the gantry baseplate, ensuring that the four (4) mounting holes align with the holes on the baseplate of the gantry.

STEP 2. Place four (4) 18-8 SS Button Head Screws, M5, 15mm through the aligned mounting and baseplate holes. See Figure 2 for reference.


Figure 2.

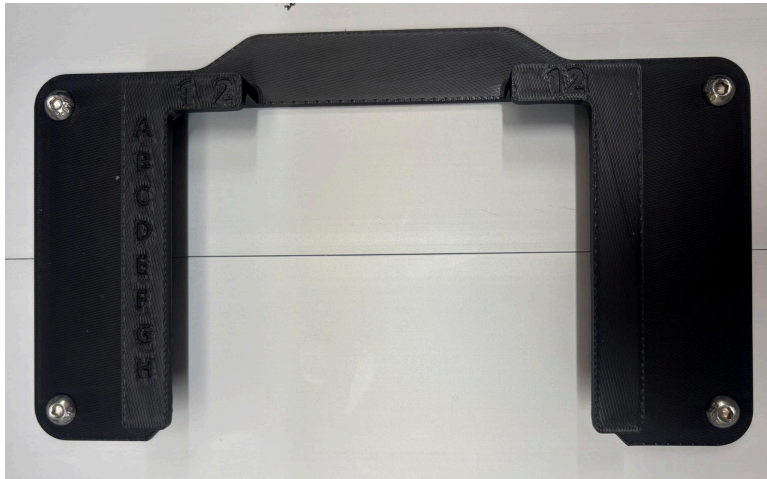

NOTE: Figure 2 depicts REV. 0 of the CubXL - ASMI Wellplate Holder.

STEP 3. Fasten four (4) Zinc Plated Steel Hex Nuts, M5 onto the four (4) 18-8 SS Button Head Screws, M5, 15mm from underneath the baseplate of the gantry. Tighten the nuts until they are snug and the mount does not move. See Figure 3 for reference.

Figure 3.

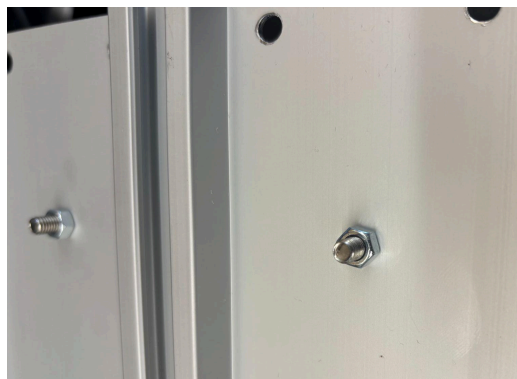

NOTE: Figure 3 depicts the underside of the gantry baseplate with the nuts installed on the screws.